

Video Processor



EOS-H-01



EA100


7.5Fr


SINGLE USE FLEXIBLE URETERORENOSCOPE

US27F-12-EU | US27F-12-US

9.5Fr ureteral access sheath

Slim without sacrificing performance of
1.2mm round working channel

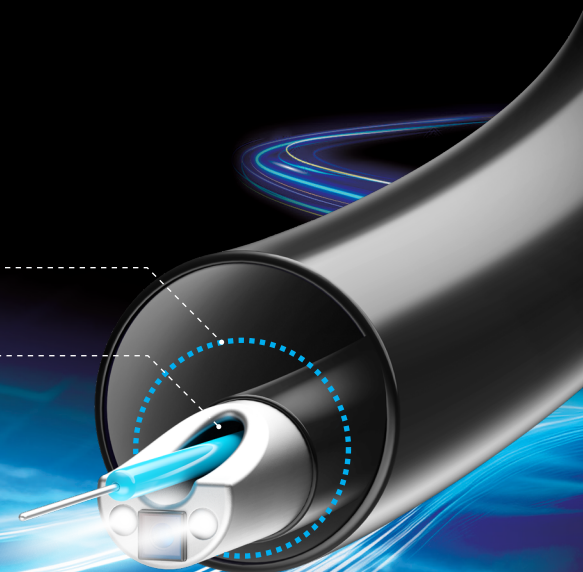

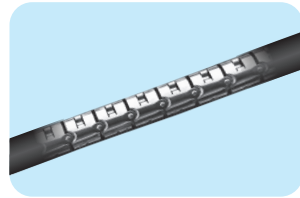

- Metal snakebone structure provide exellent maneuverability
- Protect working channel (round shape) from squeezing during deflection

Insertion section length: 670mm

- Without any adhesive tape for connection, full length metal snakebone provide consistent deflection performance during operation

275° Dual deflection



- Passive deflection function make difficult operation easy

**BRIGHT
& CLEAR**

HIGH IMAGE QUALITY



1.2mm round working channel (kink-resistant multi-layer polymer material)



Inner layer PTFE improves smoothness for passing accessories



CMOS 160K



Dual illumination



Field of view: 120°±10%
Depth of field: 3mm ~ 100mm